



Powered by **ARSOFT**

LAIA XR v.1.0.0

QUICK START GUIDE



Technical Support:
support@laiaxr.com



Legal Address &
Operational Site:
C/ del Adaja 10,
37185 Villamayor - Spain



Telephone:
(+34) 923 49 65 22



Scan this QR
code and contact
us at laiaxr.com

Edition N°: 02
Publication Date: 18/09/26

CE **FDA**
0051 Cleared

DESKTOP PHASE – PC

Step 1

>>> STARTUP & QUICK ACCESS <<<



LAIA XR

User...
Password...
Login


- 1) Turn on the PC and open the LAIA XR shortcut.
- 2) Accept the healthcare professional disclaimer.
- 3) Enter your credentials and press "Enter" or "Login".

DOWNLOAD FROM PACS

- a) Access the patient's study from PACS.
- b) Download the DICOM study. *(In some cases, a compressed file must be extracted.)*
- c) An image folder will be created. *(Remember the folder location.)*

Step 2


>>> PATIENT & STUDY <<<

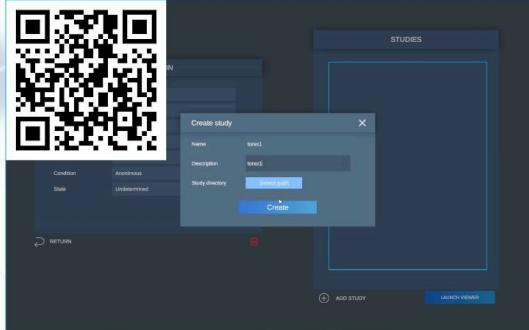

- 4) Go to "Patients" > "Add patient".
- 5) Complete the information and select the created patient.
- 6) Go to "Studies" > "Add study".
- 7) Complete the information and attach the DICOM. *(CD-ROM o USB.)*
- 8) Select the study and click "Launch Viewer".

Step 3

>>> INSIDE LAIA XR <<<

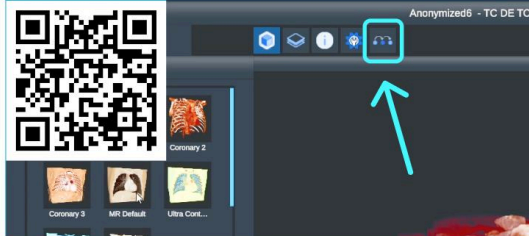

- 9) Wait for the case to load and accept the final disclaimer if it appears.
- 10) Manipulate the 3D model and slices.
- 11) Create the Surgical Planning. *(Access it using the button indicated by the arrow.)*

VIRTUAL REALITY PHASE – VR

Step 4


>>> VR SETUP <<<



- 1) Open the "Meta Horizon Link" application on the PC.
- 2) Turn on the "Meta Quest" headset.
- 3) Connect the headset to the PC using the cable and wait for the connection.

Step 5

>>> CONNECTION WITH LAIA XR <<<

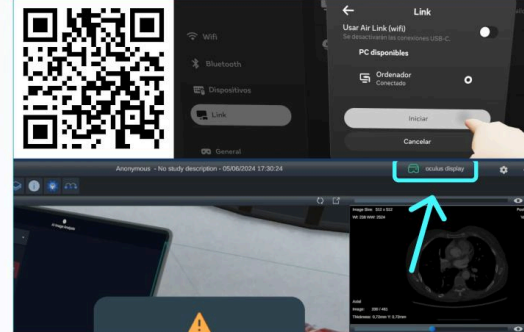

- 4) Put on the VR headset.
- 5) Connect the headset to the PC via "Link" in "Settings", inside the Meta VR environment.
- 6) In LAIA XR "Desktop", press the VR connection button.
- 7) Wait until the patient's 3D model appears. *(Do not remove the headset during the connection.)*

Step 6

>>> BASIC INTERACTION <<<



- 8) Explore the anatomy and improve your understanding through VR.
- 9) Finish the review. You can now remove the headset to exit VR mode.

SCAN THE QR CODE – ACCESS THE COMPLETE VIDEO GUIDE

Access all the Quick Guide videos and follow the step-by-step process through this QR code: from PC access to VR review.

